



**lympho-mat**<sup>®</sup>  
GRADIENT

## Multi-Jacket

### Jacket Sleeve

- Easy handling
- Comfortable to wear thanks to lightweight materials
- Can be adjusted for size using extension sets
- Can be used on a single side (12 chambers) or both sides (24 chambers)

*passion for compression*

[www.lymphamat.de](http://www.lymphamat.de)



# lympho-mat®

GRADIENT

## Multi-Jacket

### The jacket sleeve – sustainable lymphatic drainage

#### Indications:

- **Arm lymphoedema** involving the back and chest (in combination with manual lymph drainage)
- **Lipoedema**


#### Contraindications

- decompensated cardiac insufficiency
- extensive thrombophlebitis, thrombosis or suspected thrombosis
- erysipelas
- severe un stabilised hypertension
- acute soft-part-trauma of extremities
- neuropathy
- occlusive processes in the sector of lymphatic drainage
- compartment syndrome
- Acute phlegmon

In addition to manual lymph drainage and the use of compression stockings, intermittent pneumatic compression (IPC) is a method recognised worldwide for treating oedema, which is constantly being further developed by innovations from the user area.

The jacket sleeve implements a seamless intermittent compression method for arm oedema including the upper body quadrants. This allows an active compression of the entire upper body. The compression build-up begins at the fingertips and flows over the arms into the back, thorax and abdomen region.

The intermittent pressure build-up and depletion stimulates the lymphatic vasomotor system – the oedema fluid is mobilised, and its removal is sustainably promoted. Thanks to the compression gradient, a possible displacement of oedema liquid into other trunk quadrants is prevented.

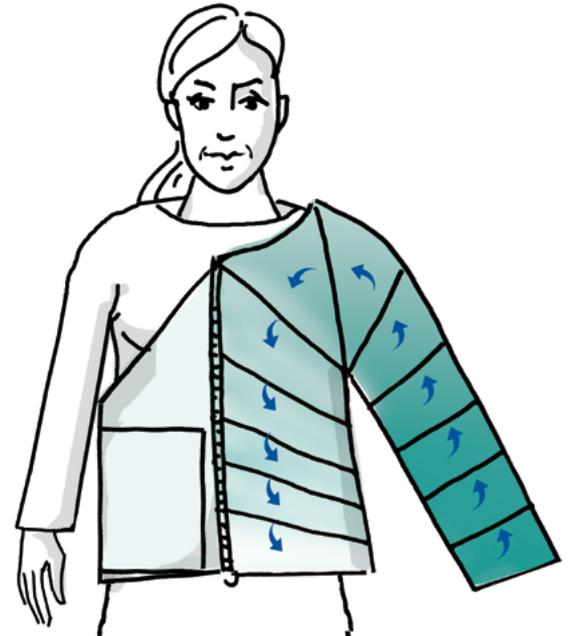

### Gradient treatment pressure

The jacket sleeve is filled with air in a compression cycle in 12 successive steps. This achieves a decreasing pressure gradient. This ensures that the mobilised liquid can disperse freely without reflux.



**12-chamber right jacket-half with bolero**  
Art.-Nr. 1182



**12-chamber left jacket-half with bolero**  
Art.-Nr. 1181



**24-chamber jacket complete**  
Art.-Nr. 1180



BÖSL Medizintechnik GmbH  
Charlottenburger Allee 13  
52068 Aachen, GERMANY  
Telefon +49 (0)241/9 00 77-0  
Telefax +49 (0)241/9 00 77-10  
info@boesl-med.de · www.boesl-med.de



**Made in Germany**